




180 liters
10,800 samples



METHER AUTOCASCADE SYSTEM ULT SERIES

-86°C ULT freezer is specially designed and manufactured for long term storage of various biological products, including viruses, germs, erythrocytes, leucocyte and cutis. Applications can be found in blood banks, hospitals, edipemic prevention service, and research institutes, laboratories in electronic and chemical plants and biological engineering institutes.

PROMED AUTO-CASCADE SYSTEM

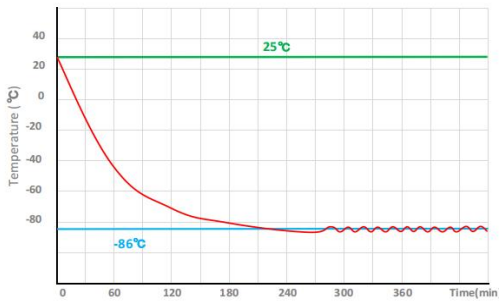

Pull down test at ambient temp.




Super-cooling technology

Auto-ascade cooling system is designed with a powerful compressor which drives the unique mixed refrigerant for fast cabinet cooling down. The system is simple in construction, energy saving and easy maintenance.



Advanced thermal insulation

Advanced insulation materials and interior design of multiple independent spaces make ULTs more efficient in Insulation and less heat loss, so as to achieve the purpose of better thermal insulation performance and energy saving.

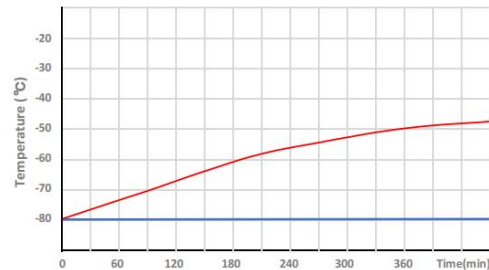

Pull up test at ambient temp. 25°C



MDF-86V838E

MDF-86V728E

MDF-86V588E

MDF-86V408E

MDF-86V340E

MDF-86V180E

MDF-86V58

58 -838 liters
2,400-60,000 samples
for options

EXTERNAL STRUCTURE

Large LED display
Full-function control panel

Ergonomic handle
Easy open
Padlock hole

Quick-release filter
Easy open & Quick clean

Proven hinges

Heavy duty castors
Universal wheels

Access port

Remote interface
RS232/485

USB interface
Power on/off

PM ULT freezers are made of color spayed anti-bacterial and anti-corrosion cold-rolled steel. Various configurations and humanized designs bring great convenience to users.

FEATURES



Large LED display
Full-function control panel



10 types of Alarms
Audible and visual alarms



Multi-layer sealing design
Better sealing, better insulation



International certifications
ISO9001, ISO13485 & CE verified




Eco-friendly
CFC-free refrigerant & insulation



Energy saving
More efficient cooling system



INTERNAL STRUCTURE

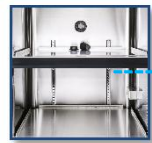


Adjustable shelves



PTC sensor



Pressure-regulated inner door handle

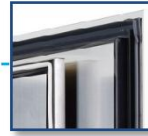

Multi-layer silicone sealing gaskets
Low temperature resistance



SUS 304 inner door filled with PURF
Detachable design

ALL
SUS 304
inside



SPECIFICATIONS

Basic		Refrigeration System		Alarms	
Cooling method	Direct Cooling	Refrigerant	R4001(mxed gas)	High/low temperature	Y
Defrost	Manual	Temp. controller	Microprocessor	Power failure	Y
Temperature range	-40℃~-86℃	System	Auto-cascade	Abnormal voltage	Y
Ambient temperature	10~32℃	Sensor	PTC 100	Sensor error	Y
Display	LED digital	Compressor brand/QTY	CUBIGEL/1	Door ajar	Y
Access port QTY/Dia.	1/25mm	Capacity(L/cu.ft)	180/6.63	Filter screen check	Y
USB Interface	Yes/1	Insulation	PURF	Thermostat failure	Y
Remote alarm Access	RS232/485	Inner door material/QTY	Foamed SUS304/3	Low battery	Y
Outside material	Color sprayed steel	Racks	9(4x3)	Condenser clean	Y
Inside material	SUS304	Boxes (2 Inches)	108	High ambient temperature	Y
Castor	Universal	Samples(2ml)	10800	Door ajar	Y
				Power failure backup (alarm)	72h
Electrical Data		Dimensions		Accessories	
Voltage frequency(V/Hz)	220V/50Hz	Package size(W*D*H)(mm)	915*928*1818	Chart recorder	/
Power(W)	660W	Exterior size(W*D*H)(mm)	830*810*1607.5	Data logger	Optional
Current(A)	3.14A	Interior size(W*D*H)(mm)	465*465*830	CO2 back up system	Optional
Power consumption	11.17Wh/24h	Shipping CBM	1.46	LN2 back up system	Optional
Noise(db)	52	N.W/G.W(kg)	207/222	Gloves	Optional

ACCESSORIES



LN₂/CO₂ backup system



Data logger



Chart recorder



Racks & Boxes



Gloves